

YAGO


▼ relax in style





TABLE OF CONTENTS

- 5** Safety instructions
- 6** Medical information
- 6** Care instruction
- 7** Warranty conditions
- 7** Start-up
- 7** Technical data
- 8** Description & operation
- 9** Remote control
- 10** Construction details
- 11** Function
- 11** Conformity declaration

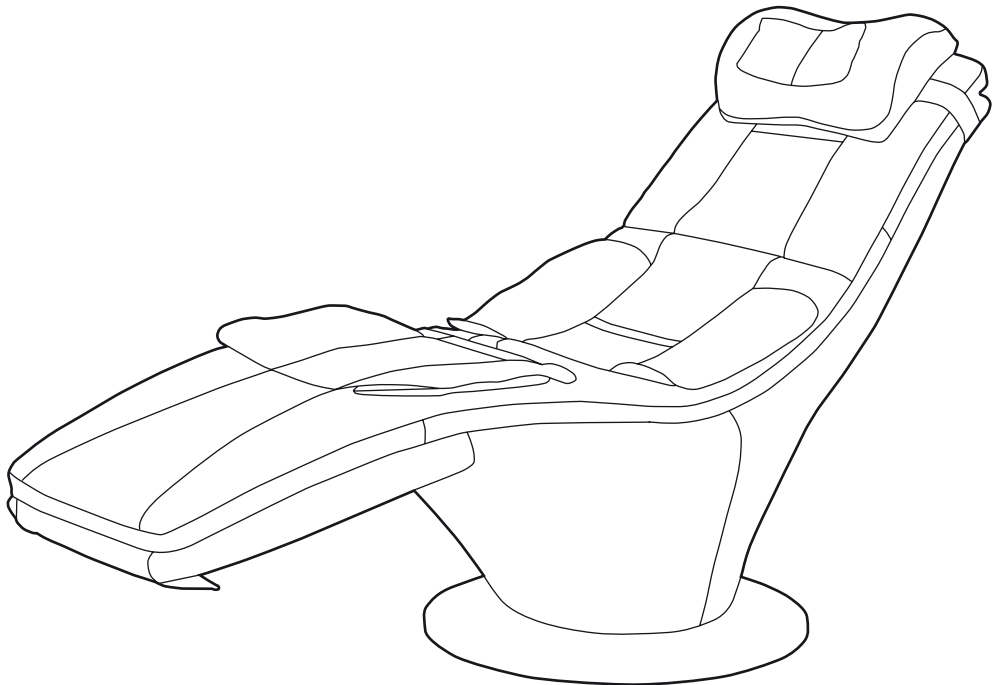



CONGRATULATIONS!

By buying your massage chair, you have proven that you care about your health. To make sure that you can enjoy the advantages of this system for a long time to come, please carefully read and observe these safety instructions.

We hope you enjoy your **YAGO**.

SAFETY INSTRUCTIONS

Read these safety instructions very carefully before using the massage chair for its intended usage. If you suffer from any of the medical conditions listed below, seek advice from a doctor before using the massage chair.

- Acute illness
- Heart disease
- Osteoporosis
- Pregnancy
- Malignant tumour
- Perception disorder
- Any condition currently under medical treatment
- Physical disability

CAUTION! To avoid accidents, make sure that you do not insert your fingers into the massage heads.

For safety reasons, you must not open up the chair yourself.

If you experience any pain when using the chair, contact a doctor immediately.

To avoid electric shocks and injuries, do not allow children to use the chair without supervision.

Do not place or use the chair in areas that are exposed to water or high temperatures.

Do not use damaged plugs or cables or loose sockets. If plugs or cables are damaged, they must be replaced by the manufacturer, a service agent, or qualified personnel.

Do not pull, twist, or kick the chair. Do not place heavy objects on the chair.

After using the chair or in the case of a long period of non-usage, make sure that the chair is not connected to the power supply to prevent electric shocks, short circuits, and fire.

If the chair malfunctions, disconnect it from the power supply immediately.

Do not touch the plug if you have wet hands.

Do not use the chair directly after eating, as doing so can cause nausea.

The duration of the massage should not exceed the recommended time of 30 minutes. This is to prevent excessive stimulation of muscles and nerves.

MEDICAL INFORMATION

As long as your doctor/orthopaedist has not expressly told you not to have back massages, you can enjoy a massage chair massage without worrying.

The strength of the massage increases if the chair is laid down. We recommend that you start off in a sitting position.

A massage should never feel uncomfortable.

Excessive force can damage the body.

Pregnant women, the infirm, and children should not use the chair. Massage can be dangerous for people with damage to the spinal column and to those recovering from an operation. If you are unsure, consult your doctor.

CARE INSTRUCTIONS

You have chosen a product manufactured using quality materials. Such a product requires appropriate care. Do not expose the chair to direct sunlight or damp.

Treat your massage chair regularly with the recommended products to ensure its long life.



WARRANTY CONDITIONS

The dealer observes a lawful warranty period for the product described in this manual. The warranty period is calculated from the date of purchase. The date of purchase is the date stated on the sales slip.

The dealer will repair and replace faulty products acquired in Germany free of charge. This does not include wearing parts such as cloth covers. The defective device must be sent back before the end of the warranty period along with the sales slip.

The warranty does not cover faults that occur as a result of external influences or repair work or changes that were not carried out by the manufacturer or an authorized dealer.

The warranty granted by the vendor is limited to the repair or replacement of the product. With regard to the warranty, the manufacturer/dealer assumes

no other liability and is not responsible for damage that occurs as a result of the non-observance of the operating manual and/or improper use of the product. If an inspection of the product by the manufacturer reveals that

the warranty claim does not relate to a fault covered by the warranty or that the warranty has expired, the customer must bear the costs of the inspection and repair.

START-UP

Before using the chair, check that the power supply voltage matches the voltage specified on the mains plug. If you have not already done so, connect the mains cable to the back of the massage chair and then insert the mains plug into the socket.

TECHNICAL DATA

Dimensions: 90 x 172 x 65 cm

Inclination angle: 130°

Weight: 35 kg

Voltage: AC ~ 220 V 50-60 Hz, 120 Watt

Operating time automatic: 15 minutes


Operating time manual: 30 minutes

Certificates:




DESCRIPTION & OPERATION


- (A) pillow
- (B) side airbag (bag / shoulder)
- (C) seat-/ tight airbag
- (D) calfrest

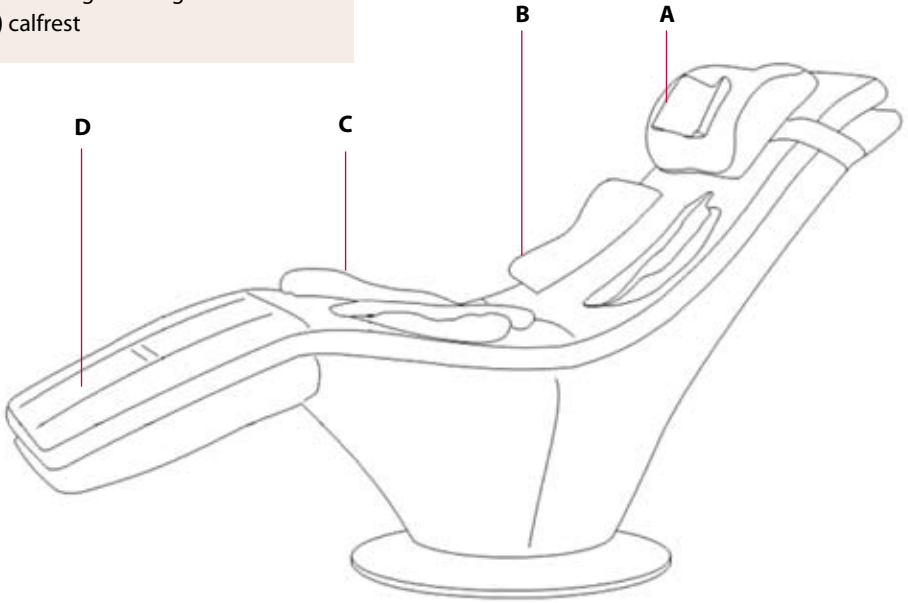


Please hold a minimum distance of 10 cm!

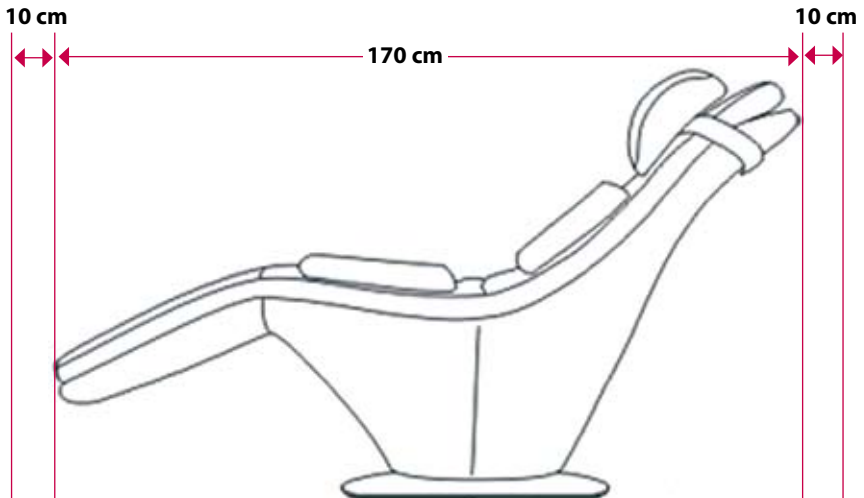

DESCRIPTION & OPERATION

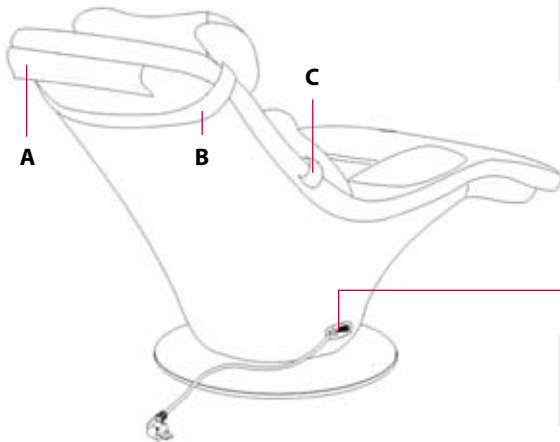

- (A)** pillow fastener
- (B)** pillow belt
- (C)** remote control holder
- (D)** energybox

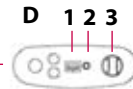




- (1)** U- Disc Input
- (2)** 3.5 mm Input for headphones
- (3)** power changer

REMOTE CONTROL







- (1)** Power on / off
- (2)** Air Massage
- (3)** Air & Vibration Massage
- (4)** Music Synchronization Massage
- (5)** Play
- (6)** Stop
- (7)** Prev / Nex

CONSTRUCTION DETAILS

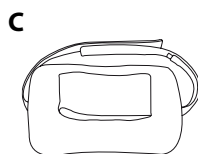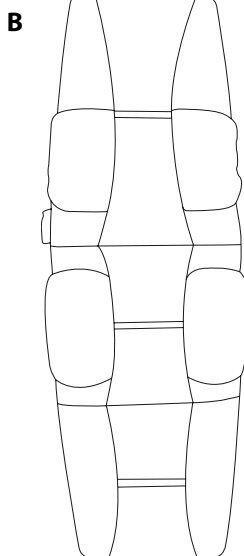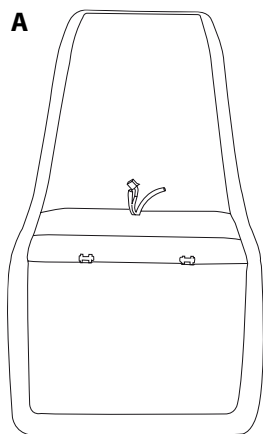

(A) Massage chair
(B) mat
(C) pillow

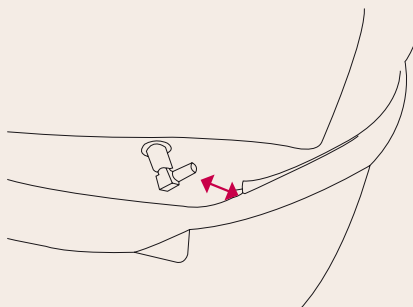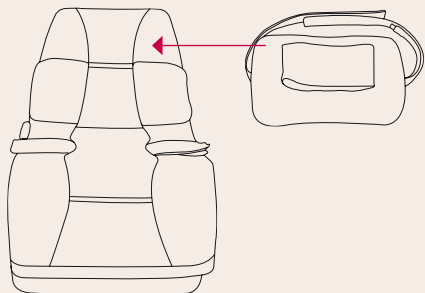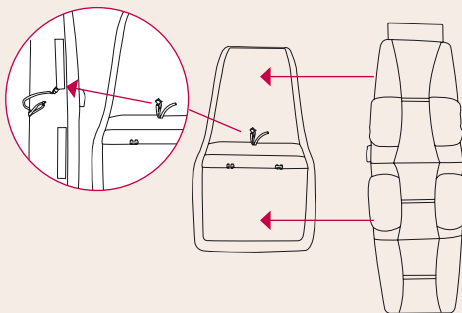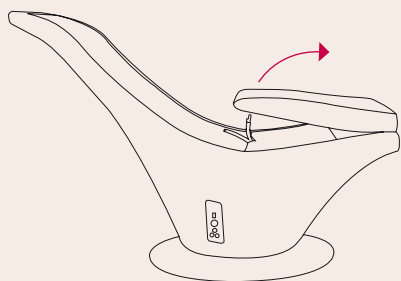

FUNCTION

High-tech meets health – The **YAGO** gives relaxation a new meaning. The unique combination of air and vibration massage is perfected with the massage synchronised to music. You determine the rhythm and speed using your own personal favourite music. Chill out and enjoy your personal massage experience every day.

CONFORMITY DECLARATION

Corresponds to the European guidelines 2006/95/EC, 89/336/EEC, 92/31EC, 93/68/EEC und 2004/108/EC.

Casada GmbH

Obermeiers Feld 3

33104 Paderborn, Germany

tel. (+49) (0 52 54) 9 30 98 - 0

fax (+49) (0 52 54) 9 30 98 - 30

Email: info@casada.de

Internet: www.casada.de

www.casada.com

Copyright

Images and texts are subject to the copyright of Casada GmbH and may not be used for other purposes without express confirmation.

Copyright (c) 2008 Casada GmbH.

All rights reserved.

time to *relax*

WITH CASADA

